

FOREO

ISSA™ 2 Range



**WORLD'S 1ST SILICONE
SONIC TOOTHBURSH**

TRAINING GUIDE



THE TOOTHBRUSH RE-ENGINEERED

FOREO is the first to apply the unique properties of silicone to oral care devices. Every device uses medical-grade silicone that is body-safe and hypoallergenic, offering safe, gentle and effective solutions.



10,000x more hygienic
than toothbrushes with nylon bristles*

*Based on clinical tests, results on file.



ULTRA-HYGIENIC

Quick-drying and nonporous to prevent bacteria buildup



SAFE

Free of phthalates & BPA



GENTLE

Incredibly smooth and soft to the touch

ELEVATE THE ART OF BRUSHING

Discover an entirely new way to brush with ISSA 2, the world's first silicone sonic toothbrush. Its sleek design and innovative features make for an exceptional brushing experience.

365 DAYS
WITH ONLY ONE CHARGE



EXTRA LONG-LASTING CHARGE

- Up to a full year of use from a single charge
- Convenient USB charging cord



VERSATILE BRUSH HEADS

- Two unique brush head designs accommodate a range of sensitivities
- Durable materials last longer than standard nylon brushes
- Fully interchangeable



SONIC PULSE TECHNOLOGY

- 2X the power
- 11,000 high-intensity pulsations per minute break up & remove plaque
- Up to 16 speeds to choose from



BRUSH BETTER WITH ISSA 2

ISSA 2 Features



Ultra-soft silicone prevents gum irritation



Leaves teeth whiter, brighter and polished



Sonic Pulse Technology clears away plaque and tartar buildup



Promotes overall oral hygiene and a healthier mouth

WHICH ISSA 2 IS RIGHT FOR ME?

ISSA™ 2



**TOUGH ON PLAQUE
GENTLE ON GUMS**

ISSA™ 2 sensitive set



2 heads are better than 1!

**THE GENTLEST
WAY TO BRUSH**

ISSA 2 RANGE



DEVICE NAME	ISSA 2	ISSA 2 sensitive set	ISSA 2 mini	ISSA mini 2 sensitive	ISSA mikro	ISSA play
# OF USES PER CHARGE	Up to 365 days	Up to 365 days	Up to 265 days	Up to 265 days	Up to 240 days	Up to 400 days (AAA battery)
BRUSH HEAD TYPE	Hybrid (Silicone + PBT Polymer)	Hybrid (Silicone + PBT Polymer, Silicone)	Hybrid (Silicone + PBT Polymer)	Silicone	Silicone	Hybrid (Silicone + PBT Polymer)
POWER	2x	2x	2x	2x	1x	1x
SONIC PULSE / MIN	11,000	11,000	9,000	9,000	6,000	9,000
# OF INTENSITIES / SPEED	16	16	1	1	2	1
USB RECHARGEABLE	✓	✓	✓	✓	✓	
SONIC PULSE TECHNOLOGY	✓	✓	✓	✓	✓	✓
TRAVEL LOCK	✓	✓				
SMILE HELPERS FUNCTION			✓	✓	✓	
ACCESSORIES	✓	✓	✓	✓	✓	
100% WATERPROOF	✓	✓	✓	✓	✓	✓
WARRANTY	✓	✓	✓	✓	✓	✓
PRICE	\$169	\$189	\$119	\$119	\$99	\$49
AVAILABLE COLORS	●●●●○	●○	●●●●●	●●	●●●●●●	●●●●●●●


HOW TO USE THE ISSA 2

Don't change your brushing routine, just your toothbrush!



APPLY TOOTHPASTE

Use your regular toothpaste and turn on your ISSA 2 by pressing the center button.

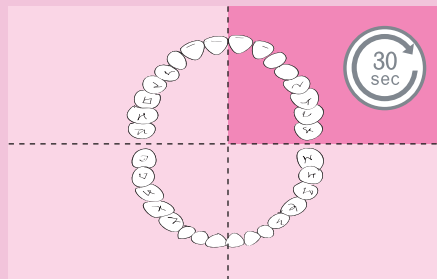

BRUSH

Brush your teeth using wide circular motions for 2 minutes; 30 seconds in each quadrant of the mouth.



RINSE

Turn off your device by pressing the same center button. Rinse your mouth and your ISSA 2 with water. Repeat twice per day for a fresher, brighter smile!

TONGUE CLEANER (OPTIONAL ATTACHMENT)


For a full-mouth clean, replace your brush head with the tongue cleaner, and use vertical strokes to gently scrub the tongue.

REPLACING YOUR ISSA 2 BRUSH HEAD

Lift the brush head straight up to detach. With the bristles of the replacement head facing forward, align its center hole with the metal dowel and push down to attach securely.

We recommend changing your ISSA 2 brush head every 6 months. Each *sensitive* brush head and Tongue Cleaner lasts up to 1 year.

All ISSA brush heads are compatible with the ISSA 2 range.

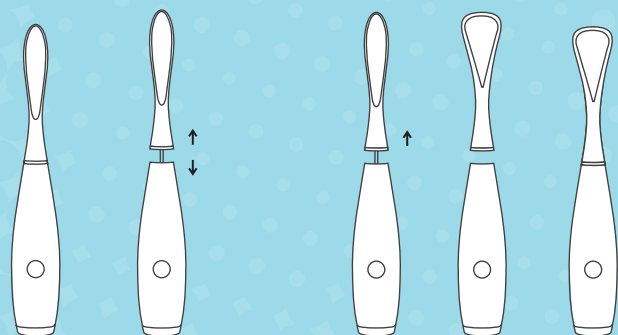

CHARGING YOUR ISSA 2

When the indicator light begins to blink, it's time to charge your device. Use the USB charging cord to connect to any outlet for 1 hour. When the indicator light stops blinking, your ISSA 2 is fully charged.

365 DAYS
WITH ONLY ONE CHARGE

COMPLEMENTARY PRODUCTS

BRUSH HEADS

Hybrid



sensitive



Tongue Cleaner



TOOTHPASTE



Try each individually, or experience the benefits of all 3 with the ISSA Professional Set

FAQs

WHAT KIND OF TOOTHPASTE CAN I USE WITH AN ISSA 2?

ISSA 2 brushes can be used with any standard toothpaste, tooth powder, tooth polish or other tooth-cleaning products. Do not use any toothpaste that contains bleach with the device. For best results, we recommend using FOREO *naturals* Smile-Perfecting Toothpaste.

FROM WHAT AGES CAN ISSA 2 BE USED?

The ISSA 2 is suitable for ages 5 and up. However, for children over 5 years of age, we recommend starting with the ISSA mini 2. For children under 5, the ISSA mikro is recommended.

HOW DO I CLEAN MY ISSA 2?

Clean the device after each use by rinsing the brush with warm water while running your fingers through the bristles. We recommend using FOREO's Silicone Cleaning Spray on the body of your ISSA 2 and rinsing with warm water afterwards for optimal results. Allow device to air-dry. We do not recommend placing the ISSA 2 in a UV sanitation case.

CAN I USE THE ISSA 2 IF I AM UNDERGOING TREATMENT FOR ANY ORAL CONDITIONS?

ISSA 2 brushes were developed for gentle yet effective cleaning of teeth and gums, and it is ideal for users who have sensitive teeth and gums, receding gums or have undergone cosmetic dentistry. However, we recommend consulting your physician prior to using an ISSA 2 if you have medical concerns.

